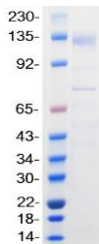


Catalog No:	LTP-V006
Product Name:	Recombinant SARS-CoV-2 Spike S1 Protein
Description:	Recombinant protein of severe acute respiratory syndrome coronavirus 2 (SARS-CoV-2) Spike S1, with a rabbit Fc tag at its C-terminus.
Alias or Clone:	2019-nCoV, COVID-19
Source:	Expressed and purified from <i>in vitro</i> cell culture of Human 293 cells
Accession No.:	NC_045512.2; YP_009724390.1; Gene ID: 43740568;
Amino acid Sequence:	The amino acid sequences of recombinant protein was derived from the Q14-A684 region of accession# YP_009724390.1
Purity:	>90% by SDS-PAGE gel and Coomassie Blue staining
SDS-PAGE & Biological Activity:	<p><i>Predicted MW of this product is ~100 kDa, however it runs higher than 130 kDa on the reduce SDS-PAGE due to a post-translational modification when expressing in mammalian cells</i></p> 
Formulation:	Protein formulated in a solution of PBS, pH7.2;
Endotoxin:	Endotoxin level is < 0.1 ng/μg of protein (<1 EU/μg)
Shipping, Storage and Stability:	The product is shipped with dry ice. Upon receipt, unopened vial can be stored at -80°C for over 12 months. Avoid repeated freeze/thaw cycles. Also the product can be aliquoted in the smaller size of working aliquots with the desired buffer and concentration, and stored at or below -20°C stable for 3 to 4 weeks.
Background:	The coronavirus Spike protein (S) is a large oligomeric transmembrane protein that mediates coronavirus entry into host cells. It contains S1 and S2 two subunits. Spike S1 mainly contains a receptor binding domain (RBD) that recognizes a variety of host cell surface receptors. S2 contains basic elements responsible for the membrane fusion. The coronavirus first binds to a receptor on the host cell surface through Spike S1 subunit, and then fuses viral and host membranes through Spike S2 subunit.

FOR RESEARCH LABORATORY TEST USE ONLY!